

**solo PAC**

# Continuous Liner Systems

Cost effective Continuous Liner systems for product contact / Isolator Transfers & waste disposal



Through Wall Cleanroom Waste Disposal Port



**SOLO** Continuous liner systems permit contained operations such as drum loading into an isolator - or drum / and or waste removal from a contaminated space with minimal loss of containment. Popular in Inert atmosphere operations – continuous liner minimises losses of the inert gas. The Solo continuous Liner systems are designed for easy use and maximum security and operator safety

**soloPAC** is available in Horizontal (Isolator) Mount as well as Vertical Product pack off or waste discharge.

Solo offer a wide range of continuous liner films from pre packed 3-layer coextruded polyethylene in 40 M to 80M pre-packs to spool loaded FDA grade film certified for Pharmaceutical product contact.

Contact the Solo Containment team +44 (0)1625 344014 for more details or email: [sales@solocontainment.com](mailto:sales@solocontainment.com)

# Solo Containment Continuous Liners are available in a range of sizes and formats to suit most applications :




**solo PAC**

- 40M & 80 M long Pre-Pack liners for Isolator bag-out and waste applications (Anti-Static & food grade)
- 10M & 15M Spool mounted Pre-Packs of FDA grade – anti static liner rated for Pharmaceutical product contact.

Size Range	Economy LDPE film (not Anti-Static)	FDA compliant Anti-Static film	40M / 80M pre-pack
100mm Dia.	✓	✓	
150mm Dia.	✓	✓	
200mm Dia.	✓	✓	
250mm Dia.	✓	✓	
300mm Dia.	✓	✓	✓
350mm Dia.	✓	✓	
400mm Dia.	✓	✓	✓
450mm Dia.	✓	✓	
500mm Dia.	✓	✓	
600mm Dia.	✓	✓	



FDA grade Film for ATEX environments

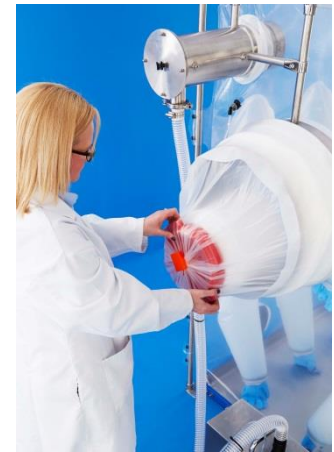

40M/80M Isolator Exit cassette

Economy / FDA quality Spool mounted films are available with high quality stainless steel mounting canisters complete with multi groove safe change over systems and quality Silicone O-rings for tension control.


Pre-pack films may be mounted on a thermoplastic carrier spool for Horizontal Mount laboratory or isolator applications . For Vertical Mount applications for dust collector waste recovery Pre-pack films are mounted on a robust stainless or epoxy coated carbon steel spool.

Custom designs are available:

### TYPICAL soloPAC LDPE film specifications

- Complies with FDA 21 CFR
- Fully compliant to 2002/72/EC and amendments
- Meets test parameters of EP 3.1.3
- Guaranteed level of antistatic.
- No animal derived components (BSE/TSE).
- Multilayer technology (PE). On pre-pack films

A full set our crimping utilities are available from Solo Containment:

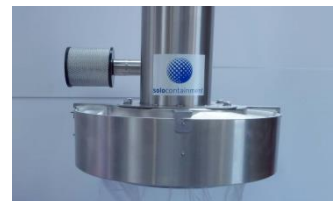

Anti-Static Dust Collector Waste Disposal



Orange Crimp System

**For expert advice contact Solo today:**

Solo Containment Ltd. Unit 6 Rupert Park, Poynton, Cheshire, United Kingdom, SK12-1PQ

www.solocontainment.com Tel +44 (0)1625 344014 Email [sales@solocontainment.com](mailto:sales@solocontainment.com)